



ARIZONA DEPARTMENT OF PUBLIC SAFETY
INTOXILYZER SERVICE REQUEST
SEND ALL COPIES TO DPS CRIME LABORATORY

QUALITY ASSURANCE SPECIALIST

| | | |
|---|---|-------------------------------|
| DATE 10/17/2008 | SERIAL NUMBER (List simulator S / N also) INTOX: 20-00003721858 SIMULATOR: | FIELD LOCATION |
| DPS EQUIPMENT NUMBER | SUBMITTING AGENCY'S FULL NAME (NO INITIALS) Goodyear PD | |
| QUALITY ASSURANCE SPECIALIST W. Will | | PHONE NUMBER (623) 8827693 |
| QUALITY ASSURANCE SPECIALIST MAILING ADDRESS PO Box 5700 | CITY Goodyear AZ | STATE ZIP CODE 85338 |

Was a S.Q.A.P. completed prior to removal from service?
 No, Explain: COULD NOT DO IT
 Yes, Date: _____ Date out of service: _____
 No Yes Do records need to be downloaded (ACTIVE C.O.B.R.A.. ONLY)

EQUIPMENT SENT:
 INTOXILYZER POWER CORD BREATH HOSE EXTERNAL PRINTER SIMULATOR KEYBOARD MODEM
 OPERATED SIMULATOR TUBE BRACKET OTHER _____

PROBLEM DESCRIPTION: (What is the instrument doing or not doing? Include observations such as processor errors or other displayed messages.)
 UNIT WILL NOT GET TO "READY" STATE TO
 START A TEST OR PERFORM MAINTENANCE

LAB USE ONLY

CORRECTIVE ACTION
 Re load software
 verify calibration & function check

COMMENTS

| | | |
|-----------------------|--|------------------|
| TECHNICIAN'S NAME | BADGE / ID NO. TITLE 6101 Elec Tech | DATE 10/21/08 |
| CALLED FOR PICKUP | DATE | CONTACT PERSON |

INTOXILYZER 8000
Location: GOODYEAR POLICE DEPT
Serial Number: 80-001858
10/21/2008

08:57:27

Standard Lot#: 808401G
Last Changed By: W. WILL #974

OAS: J. ROGERS #1123
GOODYEAR PD

Operator: J. SLOWINSKI #6101
AZDPS

Subject: TWENTY FOUR, SAMPLE, T
DOB: 03/01/1965
Sex: M Weight: 165

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 08:58:12 |
| Diagnostic Test | Pass | 08:58:41 |
| Air Blank | 0.000 | 08:59:08 |
| 0.100 Cal Check | 0.101 | 08:59:27 |
| Air Blank | 0.000 | 08:59:56 |
| Subject Test | 0.195 | 09:00:57 |
| Air Blank | 0.000 | 09:01:29 |
| Five Minute Wait | | |
| Air Blank | 0.000 | 09:06:24 |
| Subject Test | 0.195 | 09:06:40 |
| Air Blank | 0.000 | 09:07:13 |
| 0.100 Cal Check | 0.101 | 09:07:32 |
| Air Blank | 0.000 | 09:08:01 |
| Diagnostic Test | Pass | 09:08:30 |

Successfully
Completed
Test Sequence

EXHIBIT G-2
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST W.C. Will #974 AGENCY Goodyear Police Department

DATE 10/24/2008 TIME 1901

INTOXILYZER SERIAL # 80-001858 LOCATION Goodyear Police Department

- (M) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (M) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START".
- (M) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (M) 4. Air blank completed.
- (M) 5. Calibration check completed. Test results 0. 100 AC.
- (M) 6. Air blank completed.
- (M) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number: 80-001858

10/24/2008

19:00:55


Standard Lot#: 808401G

Last Changed By: W. WILL #974

QAS: W. WILL #974

GOODYEAR PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 19:01:22 |
| Cal Check | 0.100 | 19:01:41 |
| Air Blank | 0.000 | 19:02:10 |

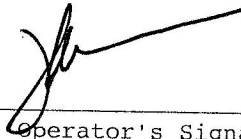

Operator's Signature

EXHIBIT G-4
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST W.C. Will #974 AGENCY Goodyear Police Department
DATE 10/24/2008 TIME 1913
INTOXILYZER SERIAL # 80-001858 LOCATION Goodyear Police Department

1. Display Reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

1. Clock time check.
 2. Date check.

OPERATIONAL TESTS

1. Alcohol-free subject test result 0. _____ AC.
 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
 4. Standard Calibration Check standard 0. 100 AC. Result 0. 100 AC.

Instrument is operating properly and accurately. YES NO _____

COMMENTS

SIGNATURE 

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/24/2008

19:11:31

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 19:13:07 |
| Diagnostic Test | Pass | 19:13:34 |
| Air Blank | 0.000 | 19:14:02 |
| 0.100 Cal Check | 0.101 | 19:14:20 |
| Air Blank | 0.000 | 19:14:50 |
| Subject Test | 0.000 | 19:15:09 |
| Air Blank | 0.000 | 19:15:36 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 19:16:54 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/24/2008

19:17:04

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | IPS* | 19:17:55 |
| Air Blank | 0.000 | 19:18:23 |

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/24/2008

19:18:32

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 19:19:21 |
| Air Blank | 0.000 | 19:19:49 |

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/24/2008

19:20:02

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 19:20:55 |
| Diagnostic Test | Pass | 19:21:22 |
| Air Blank | 0.000 | 19:21:49 |
| 0.100 Cal Check | 0.100 | 19:22:08 |
| Air Blank | 0.000 | 19:22:38 |
| Subject Test | 0.095* | 19:25:46 |
| Air Blank | 0.000 | 19:26:18 |
| Air Blank | 0.000 | 19:26:45 |
| Subject Test | XXX** | 19:27:38 |
| Air Blank | 0.000 | 19:28:12 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number: 80-001858

10/24/2008

19:29:11


Standard Lot#: 808401G

Last Changed By: W. WILL #974

QAS: W. WILL #974

GOODYEAR PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 19:29:39 |
| Cal Check | 0.100 | 19:29:58 |
| Air Blank | 0.000 | 19:30:27 |

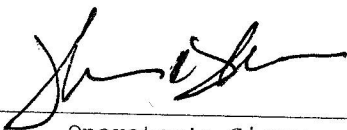

Operator's Signature

EXHIBIT G-2
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST Will AGENCY GOODYEAR PD
DATE 10/31/2008 TIME 1606
INTOXILYZER SERIAL # 80-1858 LOCATION GP

- (W) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (W) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (W) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (W) 4. Air blank completed.
- (W) 5. Calibration check completed. Test results 0. 100 AC.
- (W) 6. Air blank completed.
- (W) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE [Signature]

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number: 80-001858

10/31/2008

16:06:07

Standard Lot#: 808401G


Last Changed By: W. WILL #974

QAS: W. WILL #974

GOODYEAR PD

Test g/210L Time

Air Blank 0.000 16:06:34
Cal Check 0.100 16:06:53
Air Blank 0.000 16:07:23



Operator's Signature

EXHIBIT G-4
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST WILL AGENCY GOOSEN AR
DATE 10/31/08 TIME 1609
INTOXILYZER SERIAL # 80-1858 LOCATION SPD

() 1. Display Reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
() 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
() 2. Error recognition logic system functioning.
u Not a Successfully Completed Test Sequence printed
() 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
() 4. Standard Calibration Check standard 0. 100 AC. Result 0. _____ AC.

Instrument is operating properly and accurately. YES NO _____

COMMENTS

SIGNATURE 

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/31/2008

16:08:21

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 16:09:22 |
| Diagnostic Test | Pass | 16:09:49 |
| Air Blank | 0.000 | 16:10:17 |
| 0.100 Cal Check | 0.101 | 16:10:35 |
| Air Blank | 0.000 | 16:11:05 |
| Subject Test | 0.000 | 16:11:20 |
| Air Blank | 0.000 | 16:11:47 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 16:12:21 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/31/2008

16:12:31

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

OAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | IPS* | 16:12:57 |
| Air Blank | 0.000 | 16:13:26 |

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: GOODYEAR POLICE DEPT
Serial Number : 80-001858
10/31/2008 16:13:44

Standard Lot#: 808401G
Last Changed By: W. WILL #974

QAS: W. WILL #974
GOODYEAR PD

Operator: W. WILL #974
GOODYEAR PD

Subject: TEST,TEST,T
DOB: 11/11/1111
Sex: M Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 16:14:12 |
| Air Blank | 0.000 | 16:14:40 |

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/31/2008

16:14:46

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:15:45 |
| Diagnostic Test | Pass | 16:16:13 |
| Air Blank | 0.000 | 16:16:40 |
| 0.100 Cal Check | 0.100 | 16:16:59 |
| Air Blank | 0.000 | 16:17:28 |
| Subject Test | 0.044* | 16:20:37 |
| Air Blank | 0.000 | 16:21:08 |
| Air Blank | 0.000 | 16:21:36 |
| Subject Test | DEF** | 16:24:44 |
| Air Blank | 0.000 | 16:25:12 |

*Deficient Sample

**Deficient Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number : 80-001858

10/31/2008

16:26:29

Standard Lot#: 808401G

Last Changed By: W. WILL

#974

QAS: W. WILL

#974

GOODYEAR PD

Operator: W. WILL

#974

GOODYEAR PD

Subject: TEST,TEST,T

DOB: 11/11/1111

Sex: M

Weight: 1

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 16:27:20 |
| Diagnostic Test | Pass | 16:27:48 |
| Air Blank | 0.000 | 16:28:15 |
| 0.100 Cal Check | 0.100 | 16:28:34 |
| Air Blank | 0.000 | 16:29:03 |
| Subject Test | XXX* | 16:29:21 |
| Air Blank | 0.000 | 16:29:54 |

*Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

EXHIBIT G-2
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST Wick AGENCY GOODYEAR PD
DATE 10/31/2008 TIME 1630
INTOXILYZER SERIAL # 80-1858 LOCATION GPD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution .100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. .100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE _____

INTOXILYZER 8000

Location: GOODYEAR POLICE DEPT

Serial Number: 80-001858

10/31/2008

16:30:12


Standard Lot#: 808401G

Last Changed By: W. WILL #974

QAS: W. WILL #974

GOODYEAR PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 16:30:40 |
| Cal Check | 0.100 | 16:30:58 |
| Air Blank | 0.000 | 16:31:28 |

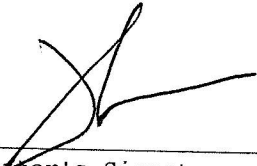

Operator's Signature